

The Vetlab SP21

Fish Disease Identification Microscope



Specialists in Veterinary
Laboratory Supplies

The Vetlab SP21 LED Monocular Microscope is the perfect, cost-effective choice for Aquarium and Koi keepers to keep their fish clean and free from parasites. Most fish parasites commonly need magnifications of x10 and x400 and the SP21's three RMS standard objectives (x4, x10, and x40), combined with the x10 and x20 eyepieces supplied, will give a magnification range of x40 to x800. A Fish Disease Diagnostic Kit is included in the purchase price and an Aluminium Carrying Case is available as an optional extra.

Features:

Optical System:

FINITE

Viewing Head:

Monocular, rotatable 360°

Eyepiece :

Wide Field (WF) eyepieces, one each of 10X/18mm and 20X/18mm

Nosepiece: Triple

Objectives:

Finite Achromatic RMS standard objectives 4X, 10X, 40X (S)

Focus Controls:

Coaxial coarse and fine focus controls. Tension adjustment on the focus mechanism. Focus Safety Stop.

Stage:

Mechanical stage with a spring clip single slide holder. Coaxial drop down controls and double Vernier scales.

Condenser:

Full sub-stage Condenser with Iris Diaphragm and rack and pinion height adjustment.

Light Source:

Built-in LED light source. Rheostat controlled.


Included: User Manual, mains power lead, spare fuse.

Optional Extras: Aluminium Carrying Case

Fish Disease Identification Kit

The Identification Kit contains microscope slides, microscope coverslips, stains, mountants and everything needed to obtain skin scrapes and gill swabs, together with a booklet of instructions and full colour photographs of most of the fish parasites that can be encountered.

For those who already have a microscope, this kit is also available as a separate item.

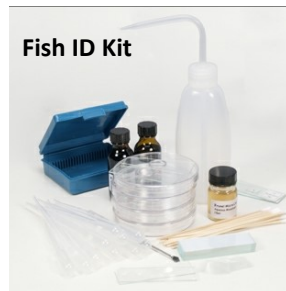

The Vetlab SP21 Monocular Microscope with Fish Disease Kit Cat. No. 1960

**Cat. No. 1954 Fish Kit only
Cat. No. 1956 Carrying Case**

2 Year Warranty

CE Certified

Visit our website to see our full range of Veterinary Equipment, Diagnostic Test Kits & Consumables...

- Special Offers
- New Products

www.vetlabsupplies.co.uk